



DICOM Export Protocol SIRONA CBCT to ANPA Planning software



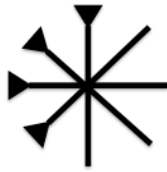



Export DICOM from SIRONA to ANPA Planning software.

1. Start Sidexis application.
2. Open the patient file and select the 3D image in Sidexis and open it (No need to click on 3D button).
3. In Sidexis software, click on **IMAGE** (on the top bar of the Sidexis Software, next to **EXAM**).
4. Select **EXPORT** in the drop down menu and click on **VOLUME FOR THIRD PARTY (DICOM)**.

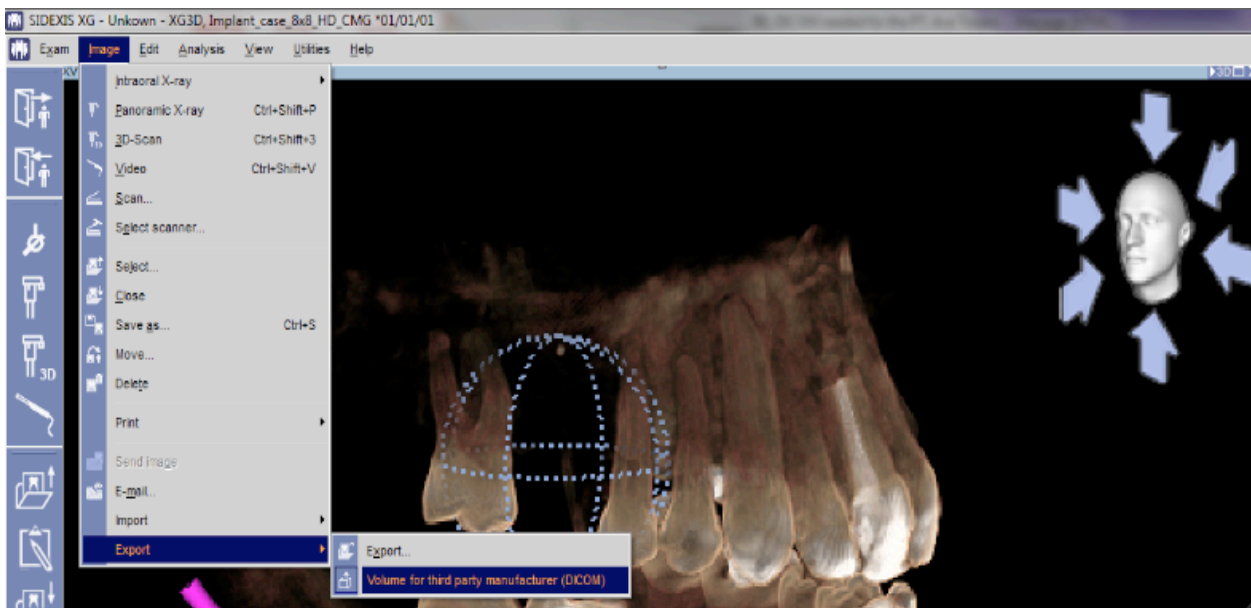


5. Select the Image area from horizontal and vertical eg: Maxiall or Mandible, and the click **OK**.

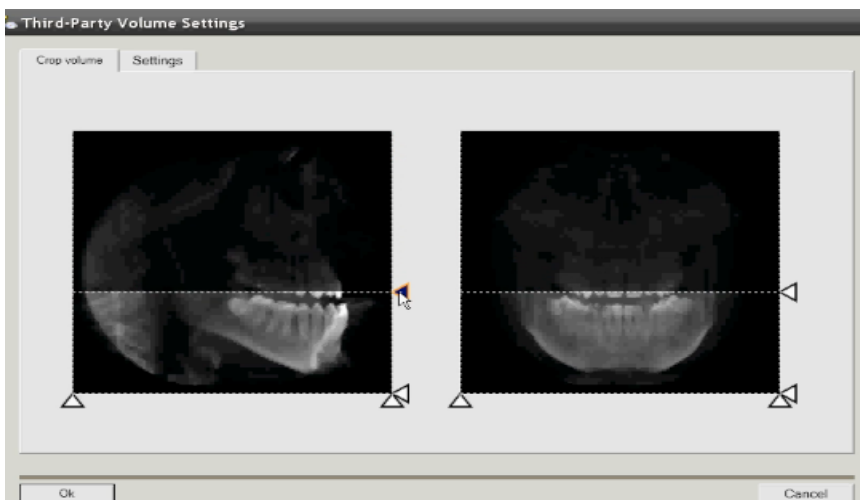


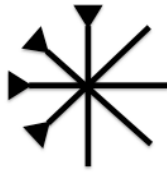


1. Select the resolution for 0.3x0.3x0.3 and click **OK**

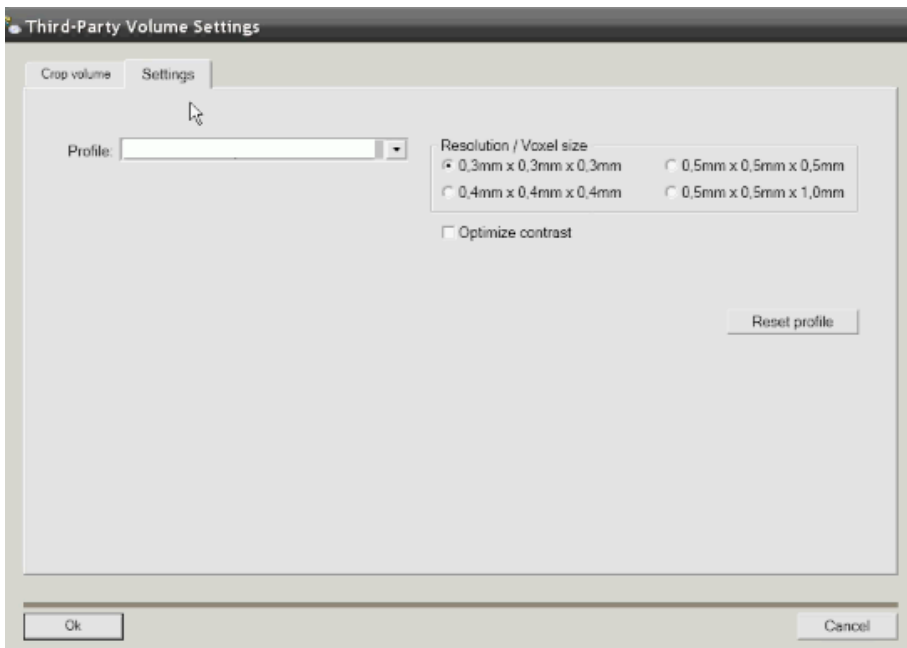

2. Select “**DICMDIR**” and choose the Export Directory “ C:\AccuGuide_CT_Data and click **OK**.

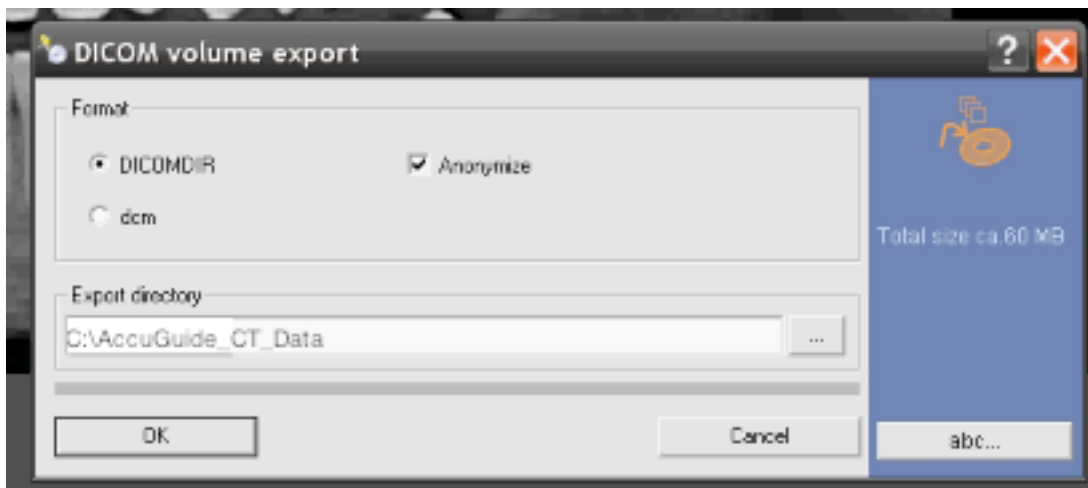

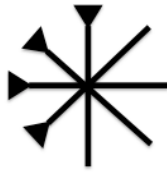

Import DICOM in ANPA Planning software

1. Start ANPA Planning software.

2. Click on  to open the DICOM.

3. Select the DICOM file and then click on Mandible or Maxilla or Head.

